



# ENDO-ONE VISION<sup>HD</sup>TRIO

▶ SDI ▶ HDMI ▶ AV RECORDER

## QUICK START GUIDE




**3<sup>in</sup>1** Video Recorder  
*for Surgical Procedures*

## SAFETY TY PRECATIONS

Before using **ENDO-ONE I visionHD TRIO**, please ensure that you read and understand the safety precautions below :

1. Do not attempt to open the case of **ENDO-ONE I visionHD TRIO** or modify it in any way, as this will void the warranty.
2. Keep **ENDO-ONE I visionHD TRIO** away from water & other liquids, as this may cause the unit to malfunction.
3. **ENDO-ONE I visionHD TRIO** may become warm during use, this does not indicate malfunction.
4. Do not place heavy objects on **ENDO-ONE I visionHD TRIO**, as this may cause the unit to overheat and malfunction.
5. If the exterior of **ENDO-ONE I visionHD TRIO** requires cleaning, wipe it with a soft dry cloth.

# Know your Remote

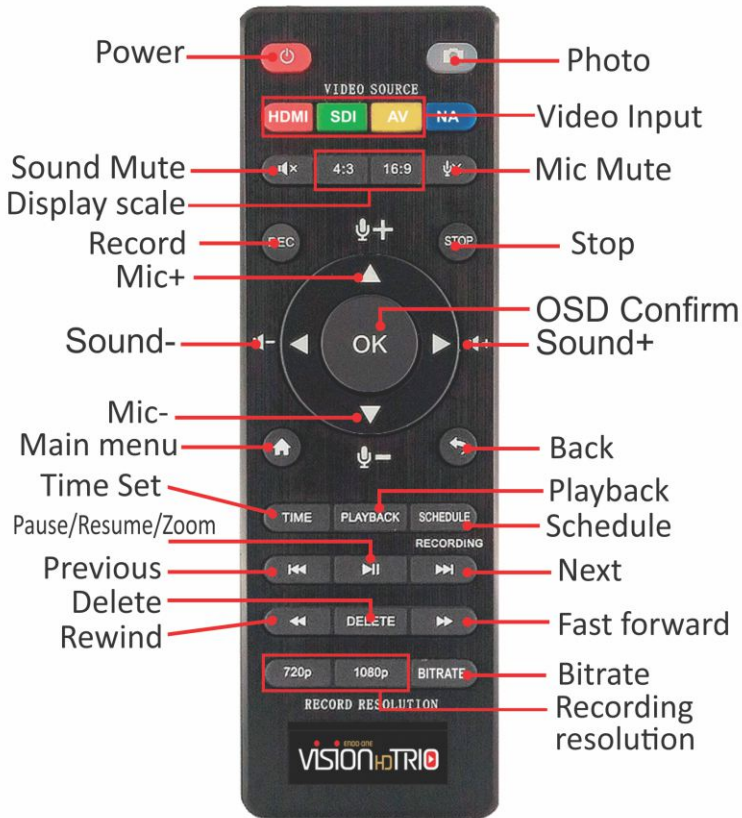

## Product Definition

No.		
1	5V/2A	Power Input
2	HDMI Input	HDMI 1.4 UP TO 4K30hz
3	AV	AV / L / R, CVBS
4	HDMI Output	Output up to 1080 60hz
5	Micro usb2.0 ( USB to PC)	Pc Mode
6	USB2.0	U disk or hard disk (32G or 2TB)
7	Line Out 3.5 mm	Headphones, speakers, etc.
8	Mic in 3.5 mm	Microphone, audio input
9	Infrared reception	
10	SDI Input	3G/HD/SD-SDI
11	USB LED	
12	SDI/HDMI/AV	SDI/HDMI/AV LED (Selection Indicator)
13	1080p/ 720p	1080p / 720p LED (Selection Indicator)
14	Power LED	
15	Resolution switching 720/1080	Selection Indicator
16	SDI/HDMI/AV Source switching	Selection Indicator
17	Record / Stop	

## External Storage Device Considerations:

1. U Disk Maximum 32GB
2. Maximum 2TB Hard Disk (If more than 1TB please format NTFS or MBR)

### **Other notes:**

1. Please use the "TIME" button for the first time to set the system time
2. After the recording is complete, wait until the light off and then remove the storage device
3. Bitrate 3 functions: Loop Recording, Limited Recording, Unlimited Recording.

**Loop Recording** : Identify U Disk or Hard Disk 12 seconds, suggest NTFS format

**Limited Recording** : Identify U Disk or hard Disk for 3 Seconds, suggest FAT32 format

**Unlimited Recording** : U Disk or Hard Disk 3 seconds, suggest NTFS format

## PC MODE:

1. The First Step is to install the software on the CD
2. The Serial Number is behind the Disc Package
3. This product can be used in addition to our software, you can also use OBS, Xsplit, etc.
4. When using the PC function, the storage device is a computer (not stored in the U disk or hard disk)





LED Light Signal Input, LED Flashing Does Not Recognize

Storage Device LED

Power

Recording LED  
(Flashing During Recording)

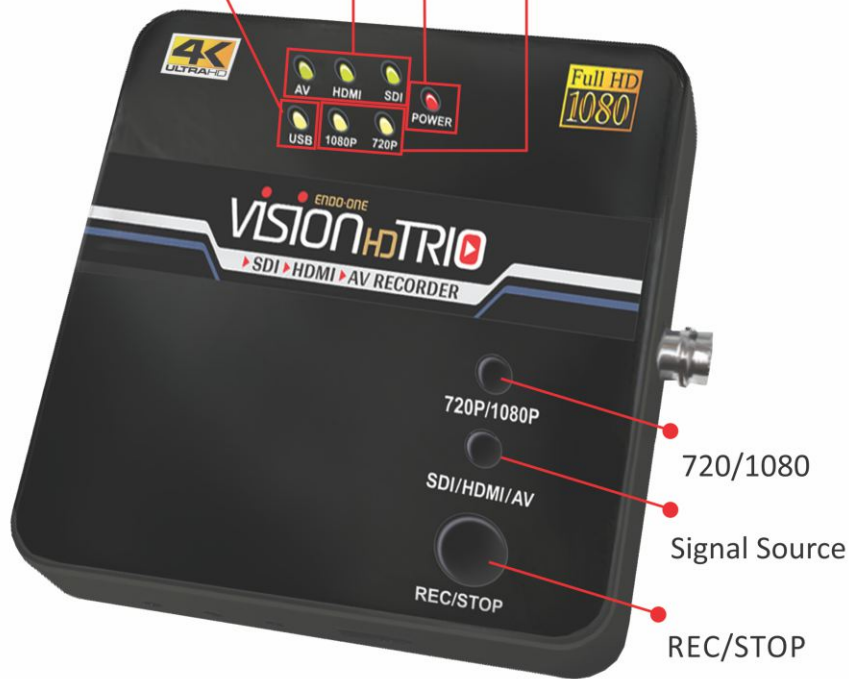

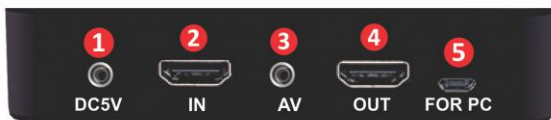





# THANK YOU

## For Choosing this Product

We hope it will provide you with years of trouble-free service.

Please look after your



If you encounter problems about using and upgrading hardware  
and software please contact our dealer or support team.

3in1

ENDO-ONE  
**VISION<sup>HD</sup>TRIO**  
▶ SDI ▶ HDMI ▶ AV RECORDER



**swift** MEDI ELEC SOLUTIONS

10/812-B, SHIVANAND NAGAR, SECTOR- 3, RAIPUR- 492008 (CHHATTISGARH)

☎ 98261 71056, 79879 19911 ✉ [swift.mes@gmail.com](mailto:swift.mes@gmail.com) 🌐 [www.swiftmedi.co.in](http://www.swiftmedi.co.in)

